


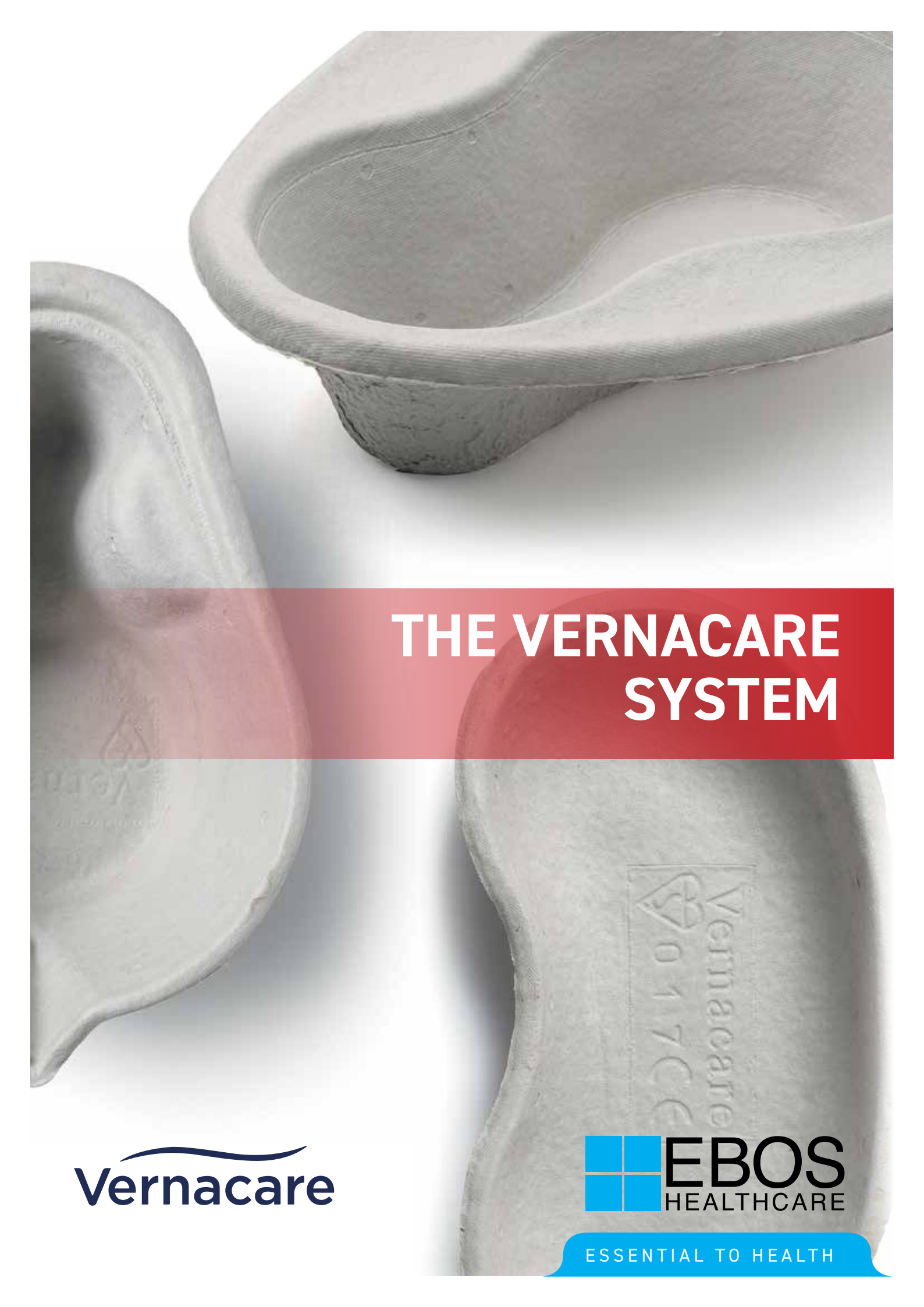

THE VERNACARE SYSTEM

**Vernacare**

**EBOS**
HEALTHCARE

ESSENTIAL TO HEALTH

Contents

Disposable Bedpan Liner and Support 4

Traditional Bedpan System 6

Urinal System.....7

Female Urinal.....8

Slipper Pan System.....9

Detergent Proof Wash Bowl 11

Bowls..... 12

General Purpose 13

Compact + Macerator..... 14

Vortex Air Macerator..... 15

Other Products you may need 16

- Tuffie 5 Cleaning and Disinfecting Wipes
- Vernagel
- Tuffie Detergent Wipes





Vernacare manufactures a complete human waste disposal system including pulp disposal unit and a full range of medical pulp products.

Experience and Innovation

Vernacare's innovative single use system and award winning products have made us global leaders in human waste management across the healthcare sector.

The system pioneered by Vernacare more than 50 years ago help reduce the risk of cross-infection and save nursing time whilst minimising the environmental and cost impact of alternative re-usable systems.

Quality Guaranteed

Vernacare's medical pulp products carry the BSI Kitemark and are manufactured in line with BS ENISO 9001: 2000 using over-issued clean newspaper. As one of the biggest pulp moulder in the world, producing around 150 million units per year. Vernacare has invested heavily in automation and quality control measures to ensure that the product is of a consistently high standard.



Vernacare

Disposable Bedpan Liner & Support



Icons Key

Capacity	
Cover available	
Stackable saves space and is easier to store	
Requires assembly	

Disposable solution for patient toileting...

Vernacare's bedpan and disposable pulp support not only helps to reduce the risk of cross infection, but also helps promote hospitals' sustainability agendas, providing environmental, financial and nursing time savings.

Products have been developed collaboratively with clinicians across the globe, to provide disposable solutions to patient toileting.

The disposable system reduces the risk of cross-infection associated with re-usable plastic products and helps tackle spiralling HCAI costs.

Bedpan Liner

Ideal for general toileting of patients and can be used either in bed or in a commode chair.

Case Quantity: 75

EBOS Code: 10021550

Supplier Code: 121AA075



Ergonomic design reduces pressure on coccyx

Curved edges allows easier patient positioning

Intuitive tear-shaped design minimises contact with patient's skin



Inner ribbed wall provides strength when used in bed

Independent support pillars distribute weight evenly

Plastic Support

To be used alongside the pulp bedpan liner with patients over 15 stone / 210lbs / 95kgs.

L400mm x W290mm x H95mm

EBOS Code: 10021195

Supplier Code: 306ZP001




Disposable Support

To be used alongside the pulp bedpan liner with patients under 15 stone / 210lbs / 95kgs.

Case Quantity: 50

EBOS Code: 10021468 Supplier Code: 122AA050


Instructions for Use:

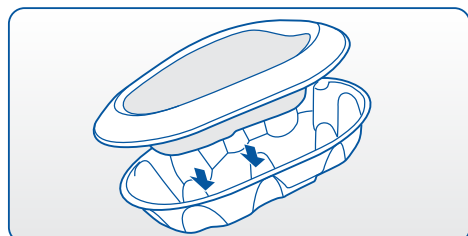


1. Place bedpan liner into disposable support, ensuring support pillars are aligned.




2. Ensure profiling bed is adjusted no higher than a 30° angle to ensure comfortable toileting.

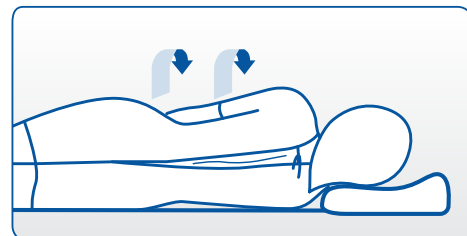

3. Gently roll patient onto their side. A second nurse may be required dependant on the patient condition. (please refer to local policy)

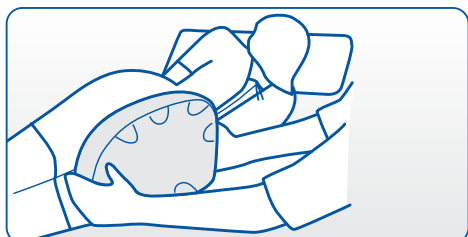

4. Position bedpan liner against patients posterior, ensuring top of the liner is aligned with patient's coccyx to ensure comfort. Gently roll patient back onto the bedpan.

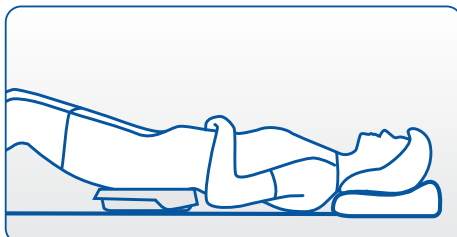

5. Ensure bedpan liner is manoeuvred into comfortable position for patient prior to toileting.



6. After toileting, secure the bedpan liner before gently rolling the patient to one side.




7. Carefully remove the bedpan, using the underside grooves to easily grip the disposable support. At this point, the disposable support can be removed and used as a lid if required.



8. Gently roll patient back into their original lying position and dispose of the bedpan liner and disposable support in the pulp disposal unit.

HOW MUCH DOES YOUR PATIENT WEIGH?

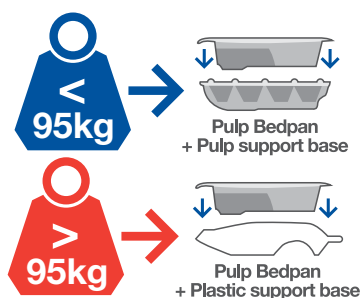

Other Products you will need:

Bedpan Liner Dispenser Rack

Wall-mounted rack removes pulp items from surfaces. Suitable for use with all bedpan liners (including NEW liners), slipper pan liners and large bowls.

L795mm x W350mm x H402mm

EBOS Code: 10020747
Supplier Code: 502ZS001



Bedpan Support Rack

Wall-mounted rack for storing up to five plastic bedpan supports. Also holds NEW plastic supports.

L625mm x W400mm x H222mm

EBOS Code: 10020746
Supplier Code: 501ZS001



Traditional Bedpan System



Bedpan Liner

Ideal for general toileting of patients and can be used either in bed or in a commode chair.

Case Quantity: 100

EBOS Code: 10021527

Supplier Code: 101AA100



Bedpan Support

Plastic support for bedpan liner, increasing strength for when used in bed.

L345mm x D305mm x H95mm

EBOS Code: 10021192

Supplier Code: 301ZP001



Bedpan Support Holder-Flat Plate

Simply clips onto the end of the bed frame for immediate access.

EBOS Code: 10020749

Supplier Code: 503ZS001



Bedpan Liner Cover

Carton of 100

EBOS Code: 10021432

Supplier Code: 197ZS100



Bedpan Support Rack

Wall-mounted rack for storing up to five plastic bedpan supports. Also holds NEW plastic supports.

L625mm x W400mm x H222mm

EBOS Code: 10020746

Supplier Code: 501ZS001



Bedpan Liner Dispenser Rack

Wall-mounted rack removes pulp items from surfaces. Suitable for use with all bedpan liners (including NEW liners), slipper pan liners and large bowls.

L795mm x W350mm x H402mm

EBOS Code: 10020747

Supplier Code: 502ZS001



Commode Attachment

Enables a disposable Bedpan to be used on a commode chair. The attachment is simply screwed to the chair in minutes.

EBOS Code: 10020754

Supplier Code: 506ZS001



Other Products you may need:

Vernagel

Vernagel is a super absorbent polymer which prevents spillages and leaks by solidifying liquids, whilst still being able to be disposed of in the Vortex pulp disposal unit.

Case Quantity: 100 6g sachets

EBOS Code: 10020733

Supplier Code: 450MA100



Electronic Fluid Measuring Scale

Designed to quickly and easily measure fluid contents inside pulp receptacle.

Automatic switch-off feature to save battery power.

Measurements in lbs and mls.

EBOS Code: 10020767

Supplier Code: 599ZS001



Urinal System



Male Urinal Bottle

Encourages independent toileting, providing convenience, dignity and comfort to patients.

Ergonomically shaped traditional design.

Case Quantity: 100

EBOS Code: 10021528

Supplier Code: 102AA100



Male Urinal Bottle (800ml) Dispenser Rack

Wall-mounted rack removes pulp items from surfaces. Suitable for use with 800ml male urinal bottles. Also compatible with female urinal.

L919mm x W400mm x H790mm

EBOS Code:10020755

Supplier Code: 508ZS001



Male Urinal Handle

Holds 1 disposable Male Urinal.

EBOS Code:10012451

Supplier Code: MUH1697



Male Urinal Bottle Holder

Rack for storing an individual male urinal at the end of a patient's bed. Suitable for use with both sizes/designs of male urinal bottle.

EBOS Code:10020752

Supplier Code: 504ZS001



Male Urinal Bottle (900ml) Dispenser Rack

Wall-mounted rack removes pulp items from surfaces. Suitable for use with 900ml male urinal bottles.

EBOS Code:10020756

Supplier Code:509ZS001



Male Urinal Lid

Suitable with all Vernacare male urinals.

EBOS Code:10021469

Supplier Code: 125AA192



Other Products you may need:

Vernagel

Vernagel is a super absorbent polymer which prevents spillages and leaks by solidifying liquids, whilst still being able to be disposed of in the Vortex pulp disposal unit.

Case Quantity: 100 6g sachets

EBOS Code:10020767

Supplier Code: 599ZS001



Instructions for Use

Use **ONE SACHET** for pulp receptacles and **TWO SACHETS** for Washbowls.

3 SIMPLE STEPS:



Female Urinal

Maximum capacity 900ml

Increased patient safety

Increased patient dignity and comfort

Case Quantity: 32

EBOS Code:10021537

Supplier Code: 108AA024



Ergonomic design helps promote self-toileting


Wide opening increases patient confidence




Curved edges allow easier patient positioning

Returned edge retains capacity and minimises risk of spillage


Instructions for Use:

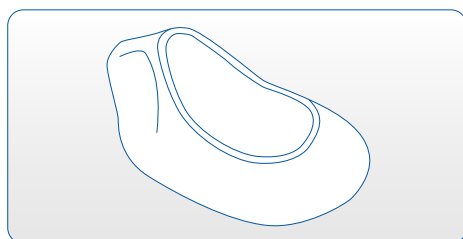


1. Product should be used with suitable patient types as appropriate.

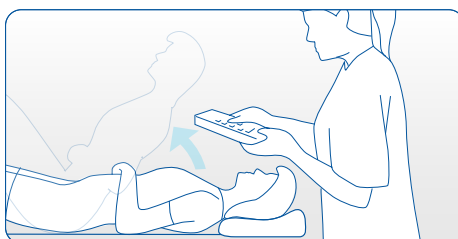


2. Adjust profiling bed to ensure patient is comfortable for toileting. Product performance will improve the more upright the patient is positioned.

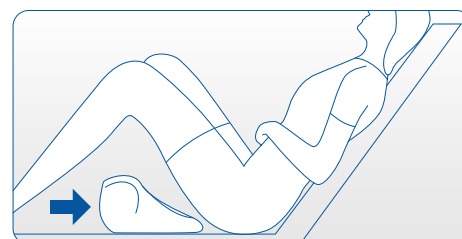


3. Position product, ensuring a comfortable fit to the patient's body. If the patient is able, they can position the product themselves; therefore increasing confidence.

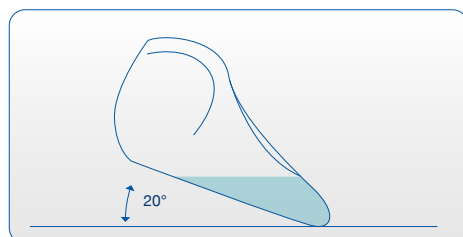

4. Product will hold up to 900ml when flat, and up to 400ml at a 20° angle.

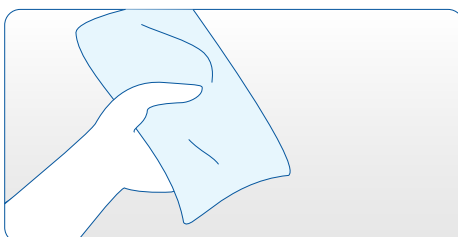

5. After toileting, secure the female urinal and use tissue paper or a soft wipe (ideally maceratable) to wipe away any excess urine.

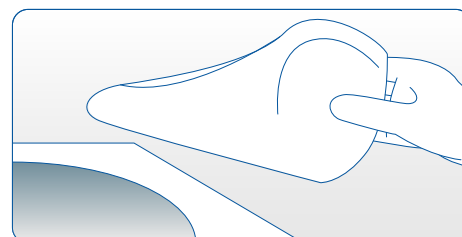

6. Carefully remove the product using the handle and dispose in the pulp disposal unit.

Electronic Fluid Measuring Scale

Designed to quickly and easily measure fluid contents inside pulp receptacle.

Automatic switch-off feature to save battery power.

Measurements in lbs and mls.

EBOS Code:10020767

Supplier Code: 599ZS001



Vernagel Bottle 475g

Vernagel is a super absorbent polymer which prevents spillages and leaks by solidifying liquids, whilst still being able to be disposed of in the Vernacare pulp disposal unit.

Case Quantity: Each

EBOS Code: 10020735

Supplier Code: 456MA475



Slipper pan System



Midi Slipper Pan Liner

Ideal for general toileting of patients with limited mobility in bed. Wedge-shaped design slips underneath patient without lifting or rolling.

Case Quantity: 100

EBOS Code: 10021538

Supplier Code: 109AA100



Midi Slipper Pan Support

Plastic support for slipper pan liner, increasing strength for when used in bed. Suitable for use with smaller slipper pan liner.

All supports are manufactured using Polyethylene.

L330mm x W285mm x H70mm

EBOS Code: 10021193

Supplier Code: 303ZP001



Maxi Slipper Pan Liner

Ideal for general toileting of patients with limited mobility in bed. Wedge-shaped design slips underneath patient without lifting or rolling.

Case Quantity: 100

EBOS Code: 10021543

Supplier Code: 114AA100



Maxi Slipper Pan Support

Plastic support for slipper pan liner, increasing strength for when used in bed. Suitable for use with larger slipper pan liner.

All supports are manufactured using Polyethylene.

L365mm x W270mm x H90mm

EBOS Code: 10021194

Supplier Code: 304ZP001



Slipper Pan Support Rack

Wall-mounted rack for storing up to five plastic slipper pan supports.

All racks are manufactured using a polished aluminium

W400mm x L625mm x H400mm

EBOS Code: 10020759

Supplier Code: 518ZS001



Slipper Pan Liner Dispenser Rack

Wall-mounted rack removes pulp items from surfaces. Suitable for use with all bedpan liners, slipper pan liners and large bowls.

L795mm x W350mm x H402mm

EBOS Code: 10020747

Supplier Code: 502ZS001



Midi Slipper Cover

Carton of 100

EBOS Code: 10021431

Supplier Code: 196ZS100



Maxi Slipper Cover

Carton of 100

EBOS Code: 10021429

Supplier Code: 194ZS100



Other Products you may need:

Vernagel Tub - 4kg

Vernagel is a super absorbent polymer which prevents spillages and leaks by solidifying liquids, whilst still being able to be disposed of in the Vortex pulp disposal unit.

Case Quantity: Each

EBOS Code: 10020734

Supplier Code: 456MA400



Electronic Fluid Measuring Scale

Designed to quickly and easily measure fluid contents inside pulp receptacle.

Automatic switch-off feature to save battery power.

Measurements in lbs and mls.

EBOS Code: 10020767

Supplier Code: 599ZS001




Slipper pan System




Instructions for Use:


The Midi and Maxi Slipper pans are designed to make toileting an easier task and the roll technique has proved to be the most effective method for patient comfort. Simply roll the patient onto the pan and roll off again after toileting.

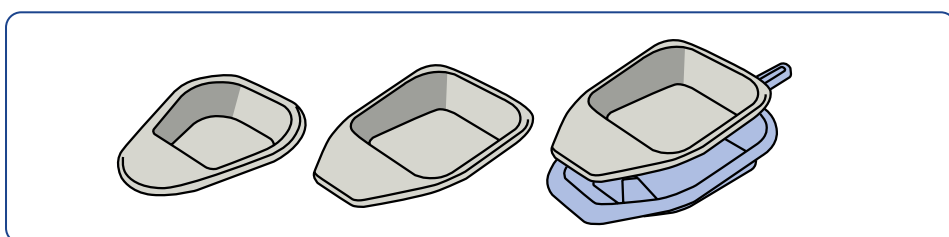


1. Choose correct size (Midi or Maxi slipper pan available) and place slipper pan liner into support.

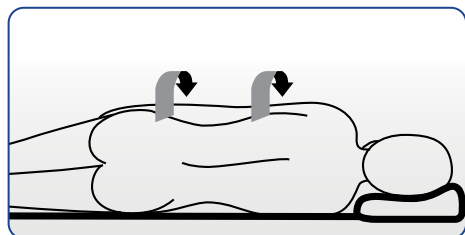


2. Roll patient onto their side.

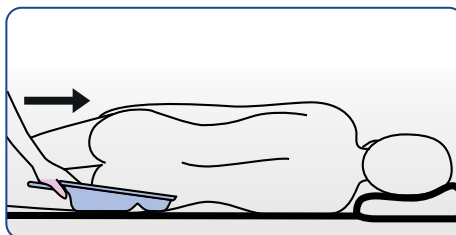

3. Carefully place slipper pan alongside patient's posterior.

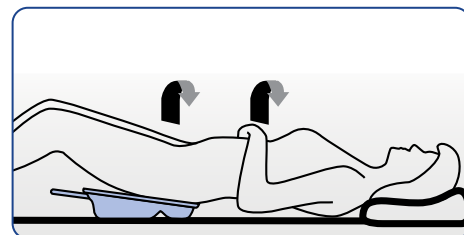

4. Gently roll the patient onto the slipperpan.

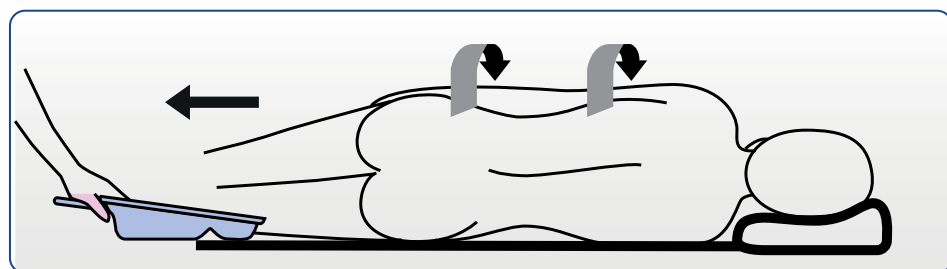


5. Roll patient onto their side to remove slipperpan and dispose of liner in the macerator.

Instructions for Use

Use **ONE SACHET** for pulp receptacles and **TWO SACHETS** for Washbowls.

3 SIMPLE STEPS:

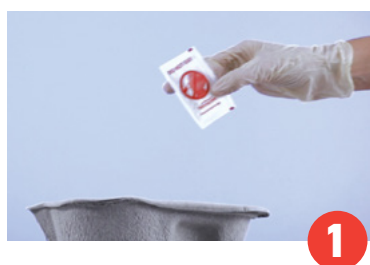

1



2



3

Detergent Proof Wash Bowl

Maximum capacity 4 litres

Increased patient safety

Holds warm water and detergent for up to 2 hours

Fill level 2 litres

Case Quantity: 100

EBOS Code: 10021546

Supplier Code: 118AA160



Instructions for Use:

Ideal for Bed baths, hand washing, general surface cleaning, wound care, shaving, continence care and foot soak. Suitable for use with detergents including soap, shampoo, conditioner and washing up liquid.



1. Take bowl from rack.

2. Add COLD water to the bowl first, this ensures the waterproof coating stays intact.

3. Then fill with WARM water. Max water temp 45°C/113°F.

4. Add detergent, soap, shampoo or conditioner and mix well.

5. For ease of carrying, use the integrated handles provided.

6. Dispose via disposal unit (maximum 2 at one time) or local waste service.

Maxi-Washbowl:

Designed for strength and easy of use. The Maxi-Washbowl holds up to 6 litres of warm water and soap for up to 2 hours (up to 45°C).



Ideal for:

- Cleaning Leg Ulcers
- Wound Care

Case Quantity: 40

EBOS Code: 10021545

Supplier Code: 117AA040

Washbowl Lid/Cover

EBOS Code: 10021486
Supplier Code: 190ZS100



Wasbowl Dispenser Rack

W383mm x D303mm x L744mm

EBOS Code: 10020757
Supplier Code: 510ZS001



Bowls



Small Bowl

General purpose bowl, ergonomically shaped with integrated graduations for measuring liquid.

Case Quantity: 200

EBOS Code: 10021530

Supplier Code: 104AA200



Small Bowl Dispenser Rack

Wall-mounted rack removes pulp items from surfaces. Suitable for use with small bowl or measuring jug.

W259mm x H800mm x D253mm

EBOS Code: 10020753

Supplier Code: 505ZS001



Commode Bowl

General purpose bowl, ideal for handling larger volumes of liquid.

Case Quantity: 200

EBOS Code:10021540

Supplier Code: 111AA200



Large Bowl

General purpose bowl, ideal for handling larger volumes of liquid.

Case Quantity: 100

EBOS Code:10021541

Supplier Code: 112AA100



Commode Pan Dispenser Rack

EBOS Code:10020758

Supplier Code: 515ZS001



Other Products you may need:

Vernagel Bottle - 475g

Vernagel is a super absorbent polymer which prevents spillages and leaks by solidifying liquids, whilst still being able to be disposed of in the Vortex pulp disposal unit.

Case Quantity: 100 6g sachets

EBOS Code: 10020735

Supplier Code: 456MA475



Electronic Fluid Measuring Scale

Designed to quickly and easily measure fluid contents inside pulp receptacle.

Automatic switch-off feature to save battery power.

Measurements in lbs and mls.

EBOS Code:10020767

Supplier Code: 599ZS001



General Purpose



Measuring Jug

General purpose measuring jug with ergonomic handle for improved grip and integrated graduations for measuring liquid.

Case Quantity: 100

EBOS Code: 10021539

Supplier Code: 110AA100



Kidney Dish

General purpose receptacle, ideal for holding a variety of medical items including dressings, wipes or surgical instruments.

Case Quantity: 260

EBOS Code: 10021531

Supplier Code: 105AA260



CS Receiver

General purpose receptacle, ideal for holding a variety of medical items including dressings, wipes or surgical instruments.

Case Quantity: 300

EBOS Code: 10021536

Supplier Code: 106AA300



Vernatray

The single use tray makes it easier for clinicians to locate the correct syringe and medication.

Case Quantity: Ctn of 60

EBOS Code: 10021547

Supplier Code: 119AA060



Jug Dispenser Rack

Wall-mounted rack removes pulp items from surfaces. Suitable for use with small bowl or measuring jug.

W259mm x H800mm x D253mm

EBOS Code: 10020753

Supplier Code: 505ZS001



Multi Cup

Ideal for taking urine/stool/sputum samples, dispensing tablets or as a measuring scoop.

Case Quantity: 200

EBOS Code: 10021544

Supplier Code: 116AA200



Specimen Collector

Ideal for collecting and measuring urine, vomit or any mid stream sampling.

Case Quantity: 100

EBOS Code: 10021542

Supplier Code: 113BA100



Multi-cup Covers

Carton of 200

EBOS Code: 10021428


Supplier Code: 193ZS200



Make the right choice

COMPACT

COMPACT+

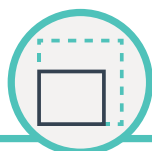

Small wards or single patient environments

Ward size & type

Large wards with multiple patient beds



Max 1 item per cycle

Items per cycle

Max 4 items per cycle



Low number of immobile patients who are toileted or washed in bed

Patient type

High number of immobile patients who are toileted or washed in bed



Patient bathroom, dirty utility room or small sluice room
Wall-mounted to save space.

Installation location

Large dirty utility room or sluice room
Compact+ is floor-standing and needs additional space for installation.



- Energy efficient
- Quiet operation
- Multiple installation options



- Minimise aerosolisation,
- spillages and splashback
- Easy to clean
- Saves space

Other Optional features include:

- Reports in real-time
- Access Control/Key
- Multiple Cycle Options (Economy/Heavy weight)

Compact
EBOS Code: 11000541
Supplier Code: 2050ZA001

Compact+
EBOS Code: 11000542
Supplier Code: 2060ZA001

Vernaclean Disinfectant and Deodoriser

VernaClean™ antibacterial deodoriser is a concentrated deodorising liquid with antibacterial properties and a fresh lemon scent.

2 Litres EBOS Code: 10020679 Supplier Code: 902SD2000



VORTEX AIR

The Vortex Air utilises the latest in motor control and diagnostics tools to minimise machine downtime whilst our SmartFlow™ maceration system ensures maximum performance. The smooth design and antimicrobial surfaces also work to reduce the risk of HCAs alongside our sensor controlled touchless technology operation.



Key Features & Benefits



Updated technology – including a deodorizer pump, capacitive sensors to measure the deodorizer and water tank levels



Smart sensors – for accurate and precise diagnostics ensuring optimal performance



Water efficient – processes four pulp items, saving on a valuable resource



Multi-cycle to suit all healthcare facilities – Eco (16L of water) Standard (23L of water) Heavy (27L of water)



Inflatable lid seal – prevents aerosolization of potentially harmful pathogens



Display – provides status of cycle during operation and time remaining to users



Multiple language options – available for international markets



Quiet operation – Only at 62.0dBA

EBOS Code: 11000437
Supplier Code: 2070ZA001

Working smarter to keep your teams safer

At Vernacare, we understand that machine downtime is not just disruptive, it can also pose an infection risk. The Vortex Air contains the most comprehensive diagnostics tools and sensors ever seen on a hospital macerator. Each unit is also tested and optimised to macerate our range of medical pulp products, ensuring your machine is running in even the most heavy-duty environments.

Why choose a single-use disposable system?



Reduces cross infection risk from reusable bedpans



More convenient and hygienic



Zero aerosols released during disposal, eradicates risk of transmission via droplets




Hands-free system reduces number of potentially infected touch points



Immediate disposal reduces potential risk of contamination spray, splash backs, or spillages during cleaning



Vernacare macerators use significantly less water than bedpan washers. In addition to this, the Vortex Air only requires cold water to process pulp products resulting in less energy consumption per product.

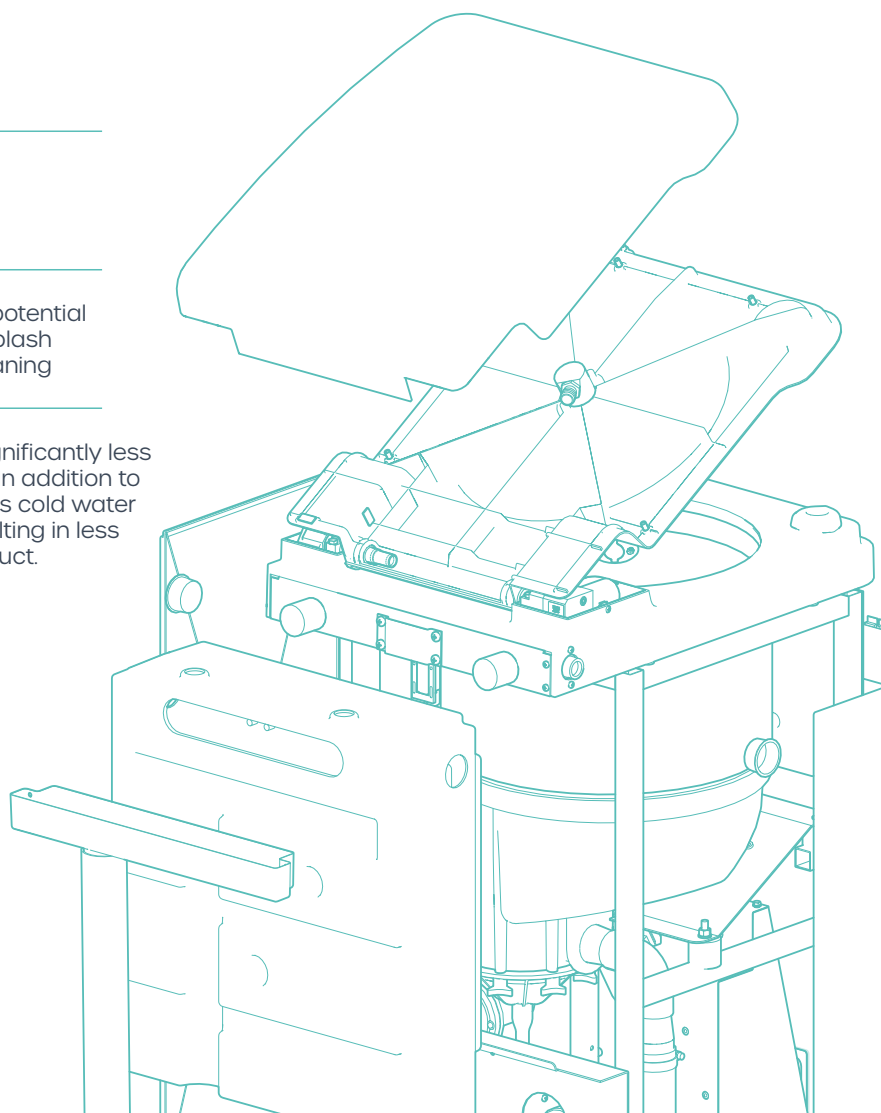

Other Products you may need:

Tuffie 5



Cleaning and Disinfecting Wipes

Tuffie 5 is impregnated with a low level biocide for cleaning and disinfection of hard surfaces and equipment.

The formulation combined with mechanical action when using the wipe aids the removal of organic matter and improves infection prevention and control.

Tuffie 5 wipes are suitable for use on all hard surfaces where cleanliness and eradication of bacteria are required.

Alcohol and Chlorine FREE!

Tub of 225 Wipes

EBOS Code:10020677 Supplier Code: 901SW225

Flexican of 150 Wipes

EBOS Code:10020675 Supplier Code: 901SW150FC

Dispenser Refill

EBOS Code: 10020676 Supplier Code: 901SW150RF

Dream Dispenser

EBOS Code:10020760 Supplier Code: 520ZS001

Vernagel



What is Vernagel?

Vernagel is a super absorbent polymer which prevents spillages and leaks by solidifying liquids, whilst still being able to be disposed of in the Vortex pulp disposal unit.

- Eases nurse workload and saves laundry costs
- Reduces risk of cross infection
- Minimises residual odours
- 1 sachet solidifies up to 1.2L of liquid.

Case Quantity: 100 x 7g sachets

EBOS Code:10020733 Supplier Code: 450MA100



Vernagel Bottle - 475g

EBOS Code: 10020735

Supplier Code: 456MA475



Vernagel Tub - 4Kg

EBOS Code: 10020734

Supplier Code: 456MA400

Tuffie Detergent Wipes

Perfect for using where general cleanliness is required quickly and conveniently.

- More convenient and hygienic than bowls of soap and water
- Can be used for cleaning equipment between patients, dressing trolleys, mattresses
- Low foaming surfactant means no need for drying cloth
- Do not dry out like alcohol wipes

Tuffie Detergent -225 Tub

EBOS Code: 10020672 Supplier Code: 901DW225

Bracket EBOS Code:10009577 Supplier Code:TTWB225

Tuffie Detergent Flexican- 6 x 150 wipes

EBOS Code: 10020670 Supplier Code: 901DW150FC

Tuffie Detergent Refill - 6 x 150 wipes

EBOS Code:10020671 Supplier Code: 901DW150RF

Dream Dispenser:

EBOS Code:10020760 Supplier Code: 520ZS001



Notes

[illegible]

Notes

[illegible]

EBOS Healthcare - Australia

300 Victoria Street,
Wetherill Park,
NSW 2164
Free Phone: 1800 269 534
Free Fax: 1800 810 257
Email: customerservice@ebosgroup.com.au
Online Ordering: www.eboshealthcare.com.au

WAREHOUSE LOCATIONS:

New South Wales:
300 Victoria Street,
Wetherill Park,
NSW 2164

Queensland:
Unit 2, 39 Distribution Street
Motorway Business Park
Larapinta
QLD 4110

Victoria:
34-78 South Park Drive,
Dandenong South,
VIC 3175

Western Australia:
2/85 Sheffield Road
Welshpool
WA 6106

South Australia:
75-83 Hardys Road
Underdale
SA 5032

Tasmania:
7A Pear Avenue
Derwent Park
TAS 7009

email: customerservice@ebosgroup.com.au
Internet: www.eboshealthcare.com.au
Customer Service Freephone: 1800 269 534

 Vernacare

 EBOS
HEALTHCARE